



# Water Quality Data at High Time and Space Resolution in the Field : Expanding Spectrophotometer Capabilities with Arduino Driven Sequential Autosampler

Bryan Maxwell, François Birgand, Kyle Aveni-Deforge and Brad Smith

## Need for both high spatial and temporal resolution

- Biological and chemical processes often exhibit high spatial and temporal heterogeneity
- Data acquisition in environmental sciences is **always partial in time and in space**
- 'Hot spots' and 'hot moments' are thought to have a disproportionate impact on biogeochemical cycles
- There is demand for tools able to **capture hot spots and hot moments together** in the field (Fig.1)
- New UV-Vis spectrometer probes make high temporal resolution water quality measurements possible for many parameters

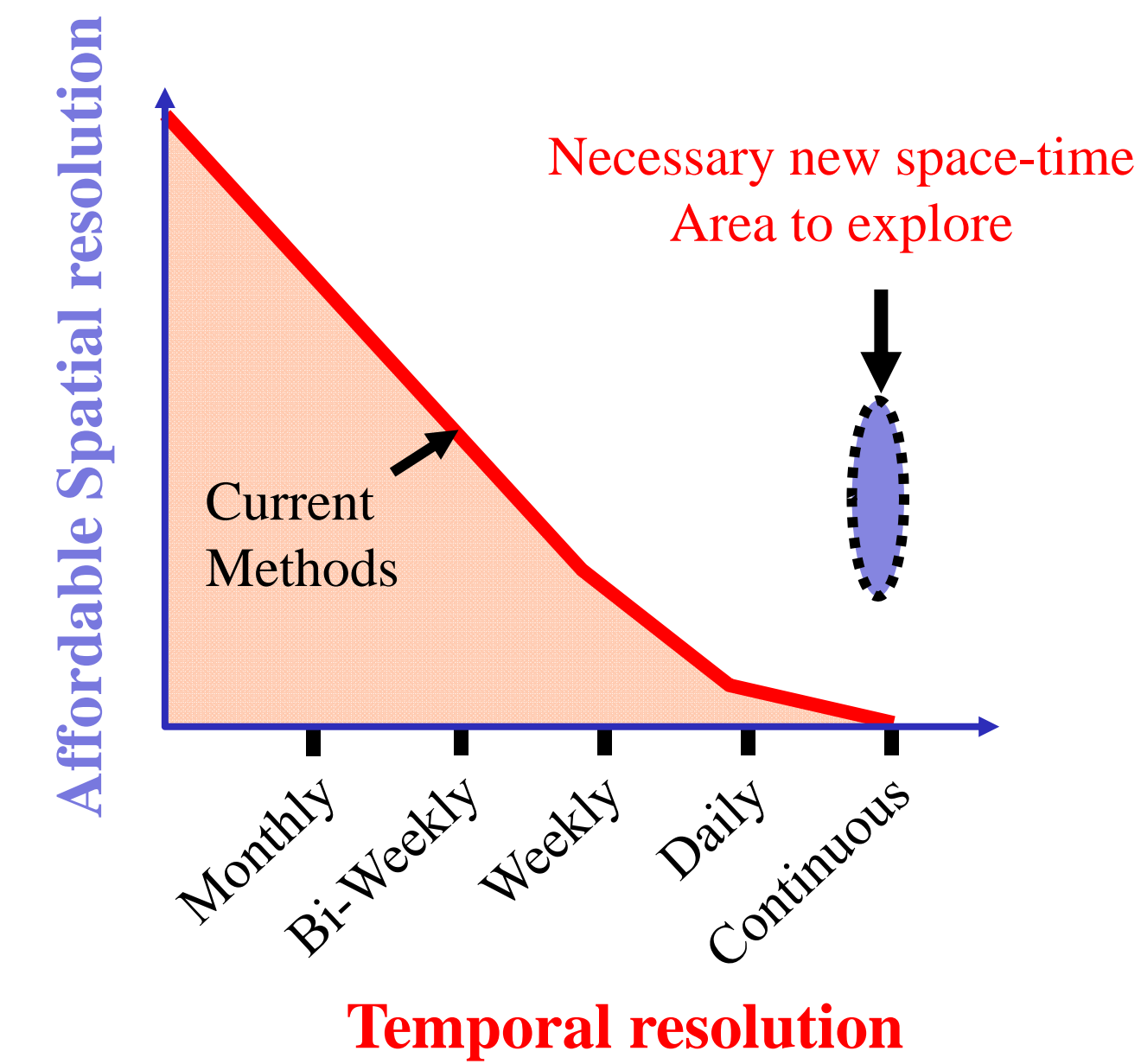


Fig. 1. Conceptual relationship between spatial and temporal resolution scales for current and desirable investigation methods in water quality and biogeochemistry

### Research Questions

- Can a single high-frequency UV-vis spectrophotometer be fitted with an autosampler to measure multiple sources of water?
- Can the volume per sample be minimized to expand use for porous media systems?
- Can such a system provide sufficient precision to address new time-space resolution demands?
- What would be the performance and applications of such a system?

## Sample Routing Configurations

- System can be configured to minimize cross-contamination
- Several internal paths are available for
  - pumping,
  - rinsing
  - purging back to source or not
  - Volume collecting
- From source to probe and back to source (A & D)
- Rinsing cycle with DI water (C & F)
- Purging previously filled line (B) can occur before sampling to probe (A)
- Purging the manifold without sending back water to source (F)
- Automatic sampling to fraction collector (E). Pump is reversed to purge the fraction collector line with air (not shown) and the left over is sent to waste (F)

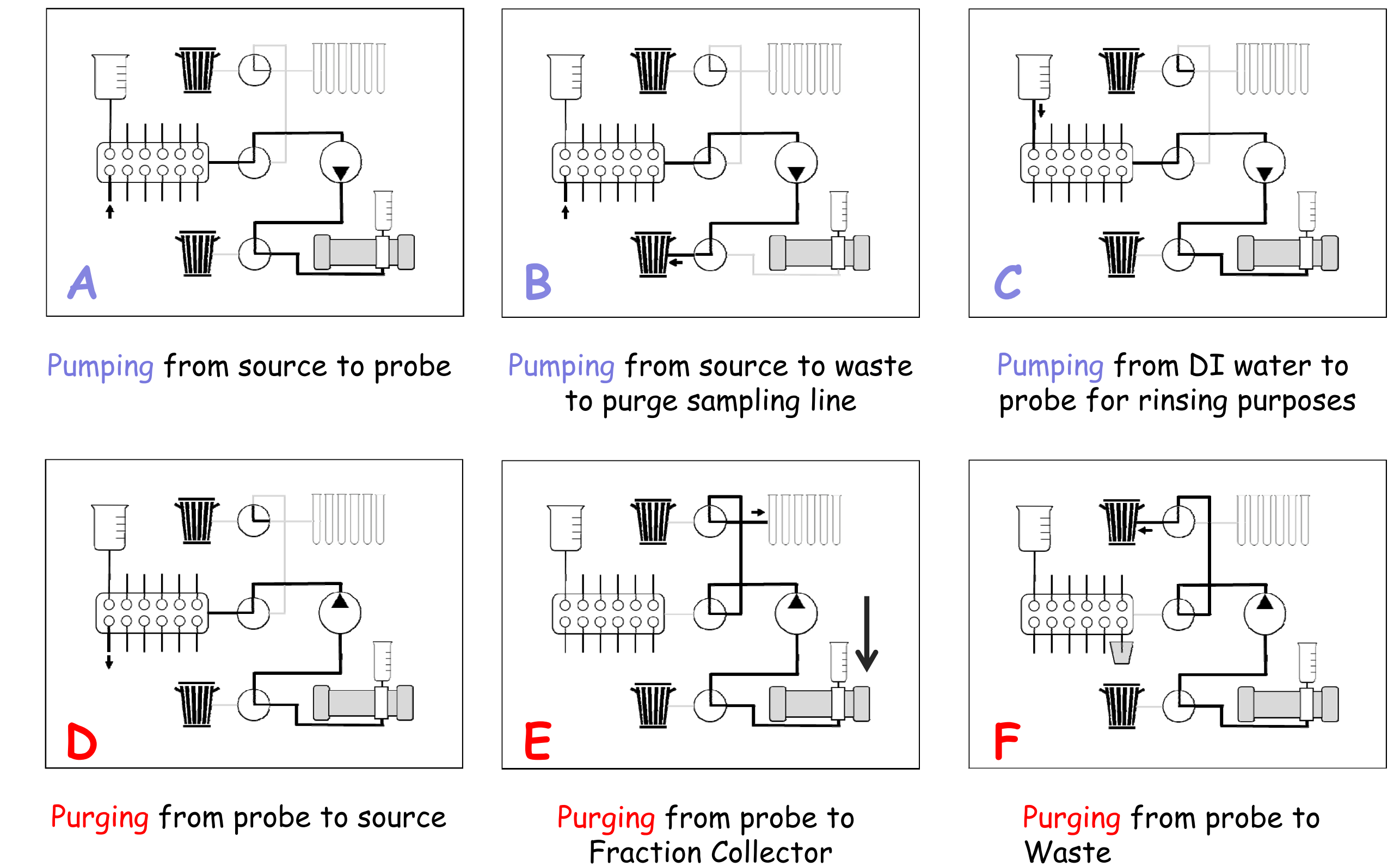


Fig. 3: Schematics of the different configurations available for pumping, purging, rinsing and sampling

## System evaluation

- Surface tension causes small water droplets to remain in 0.9 mm diameter tubing
- Potential for cross contamination
- Test system for riskiest configuration sequence, i.e. in Fig.3 **ADCF, ADCF**, etc. alternating sampling from two, low and high concentrations, bottles of 500 mL each
- Results show cross-contamination of less than 80  $\mu$ L between consecutive samples (Fig. 5)
- It would take **2,500 samples to dilute 20L** of water with DI water from 5.00 to 4.95 mg N/L of nitrate, i.e. **~ 7 days**

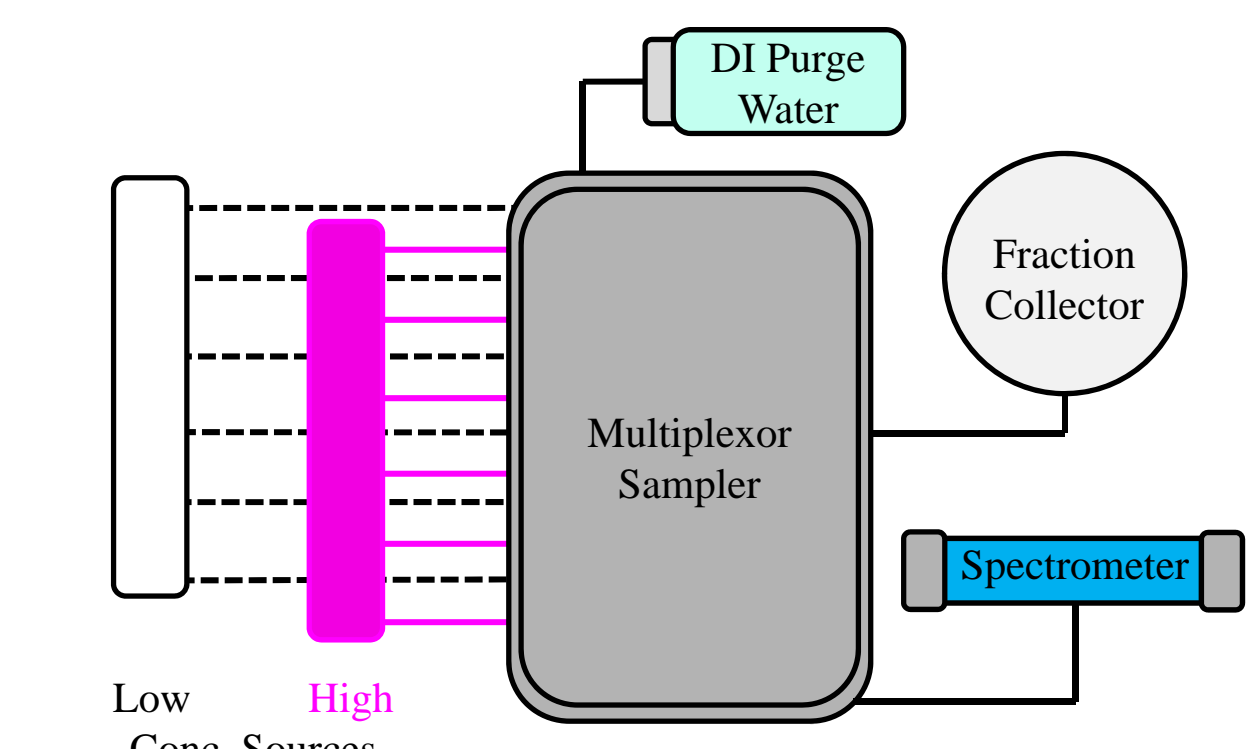


Fig.4: Alternating between low and high concentration sources to evaluate cross-contamination



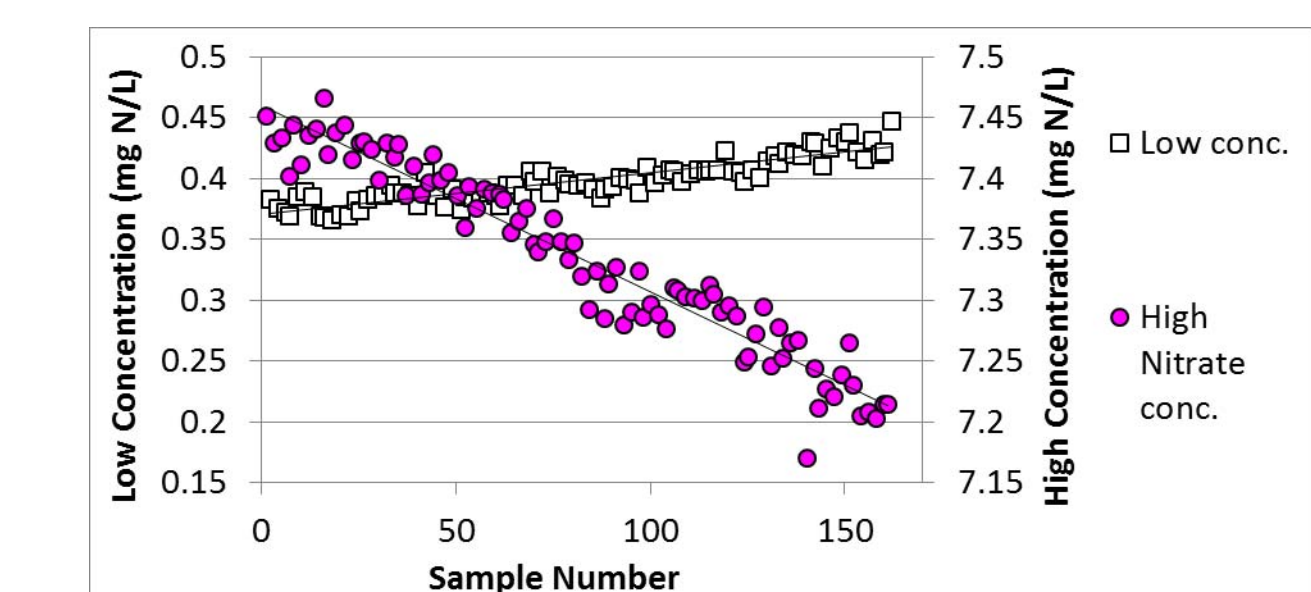


Fig.5: Changes in concentration in the two sources over time allows calculating apparent contaminating volume

## System Design and Functioning

- Small sample volume:** less than 15 ml thanks to low volume cuvette and 0.9 mm ID PTFE tubing
- Multi-port intake:** 12-way valve allows 12 different intakes from separate sources up to 12 m away
- Several automated configurations** for pumping, purging and rinsing
- Time resolution** per source better than 1 hour for each of the 12 sources
- Arduino controlled** fully automated; works on 12V DC
- Affordable sampler:** less than \$5,000 to build
- UV-vis spectrometer** : gives immediate, and rich data stream which can be analyzed to provides accurate results at high sampling frequency.
- Fraction volume collector** : collects samples for lab analysis to calibrate/validate spectrometer results

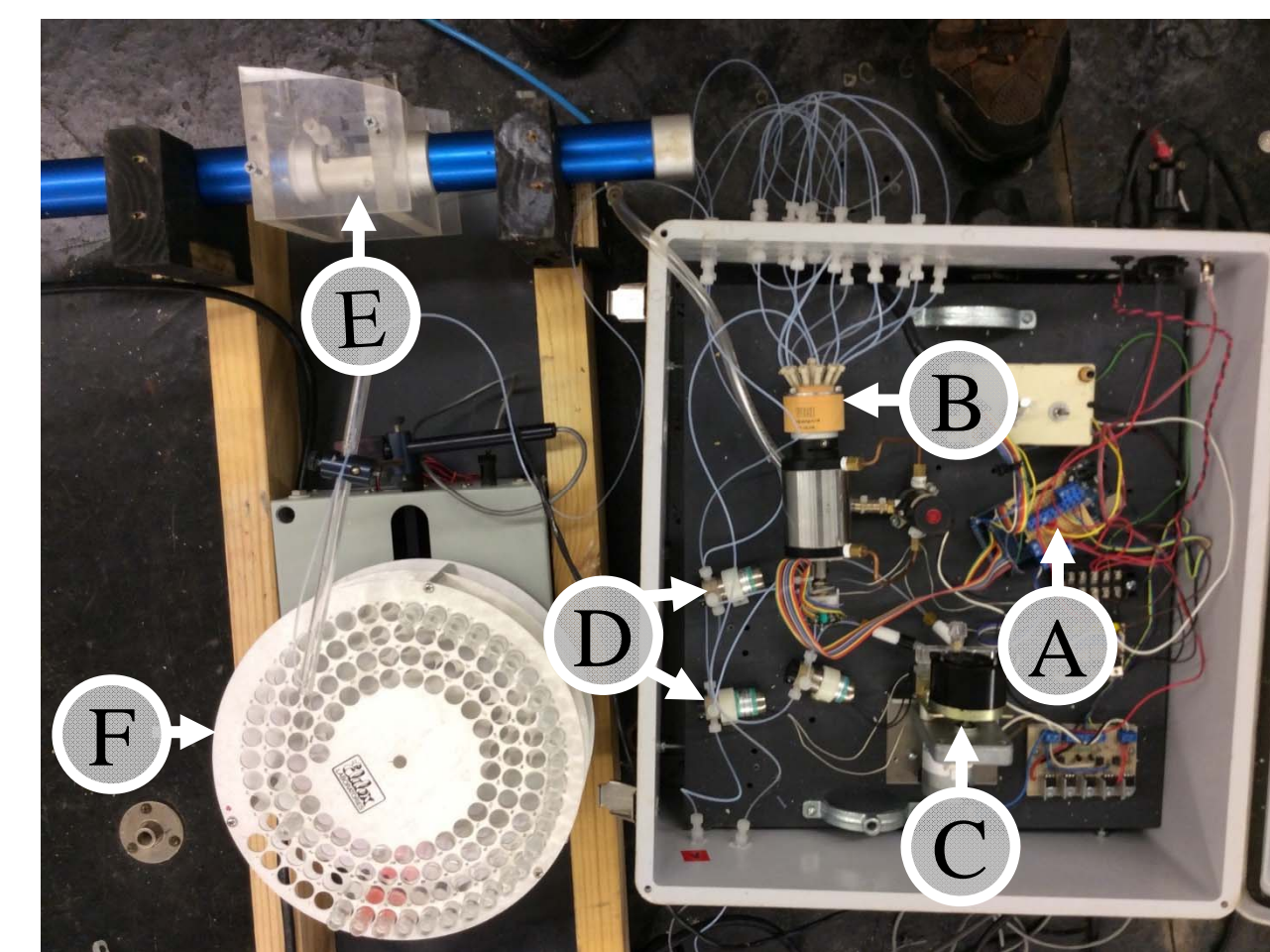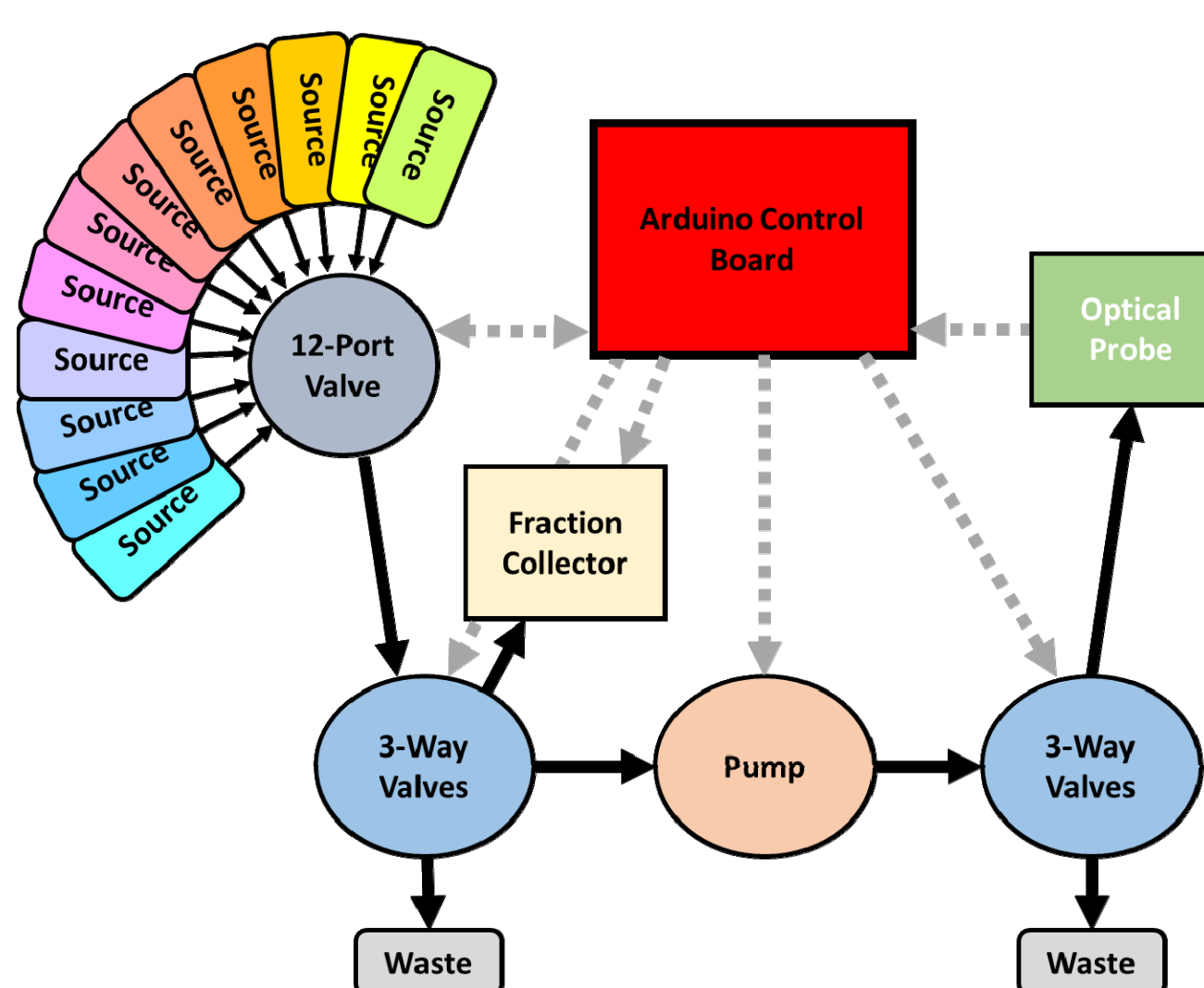



Video of the system At work




Fig. 2. Top: flow chart of system. Bottom: An Arduino control chip (A) controls selection of the 12-way valve port (B), the pump (C) and 3-way valves (D) to control sampling configuration to spectrometer (E) to fractional volume collector (F) for lab comparison.

## Application: in stream mesocosm study with replication

- 6 min per sample
- 36 min resolution for each mesocosm
- Allows refined kinetics calculations

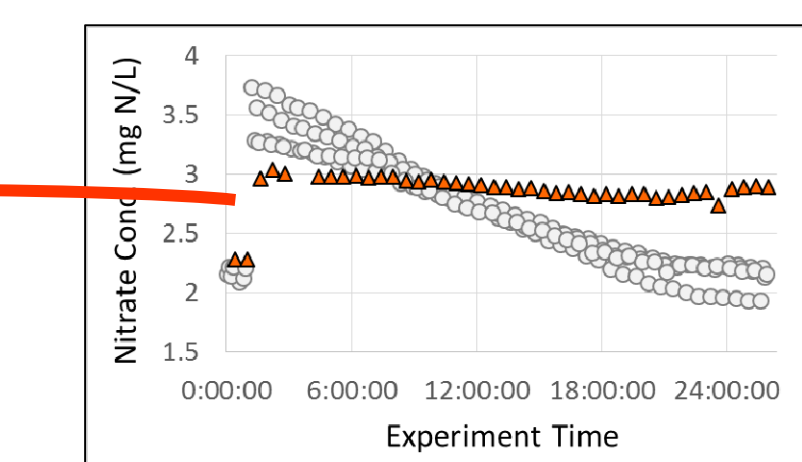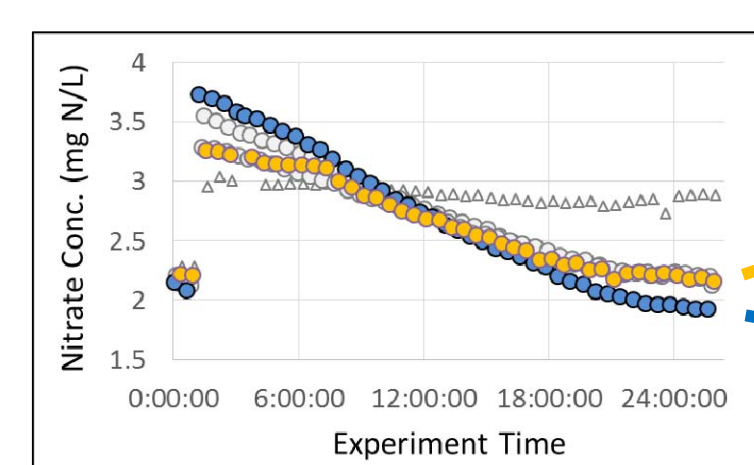

Control

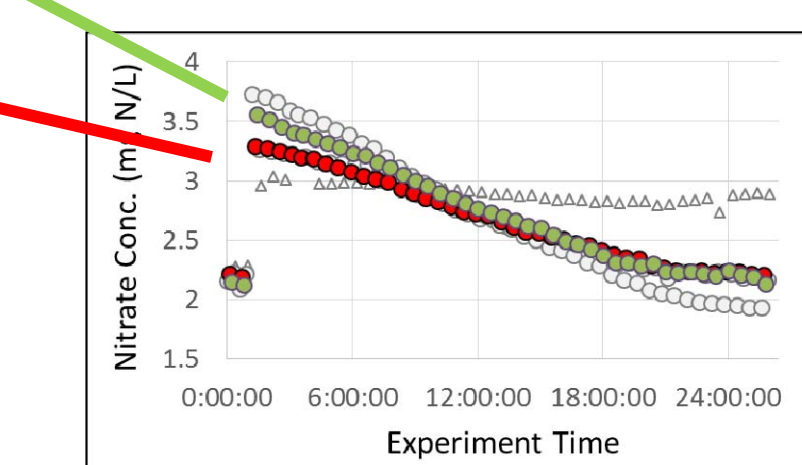

## Conclusion

- Mini-volume automatic analyzer opens the possibility to explore new time-space resolution domain
- Low Cross-contamination warrants use of instrument in most situations
- Promising method for studying 'hot spots' and 'hot moments'

Please visit our lab:

